

DENTAL UNIT CORRECTA








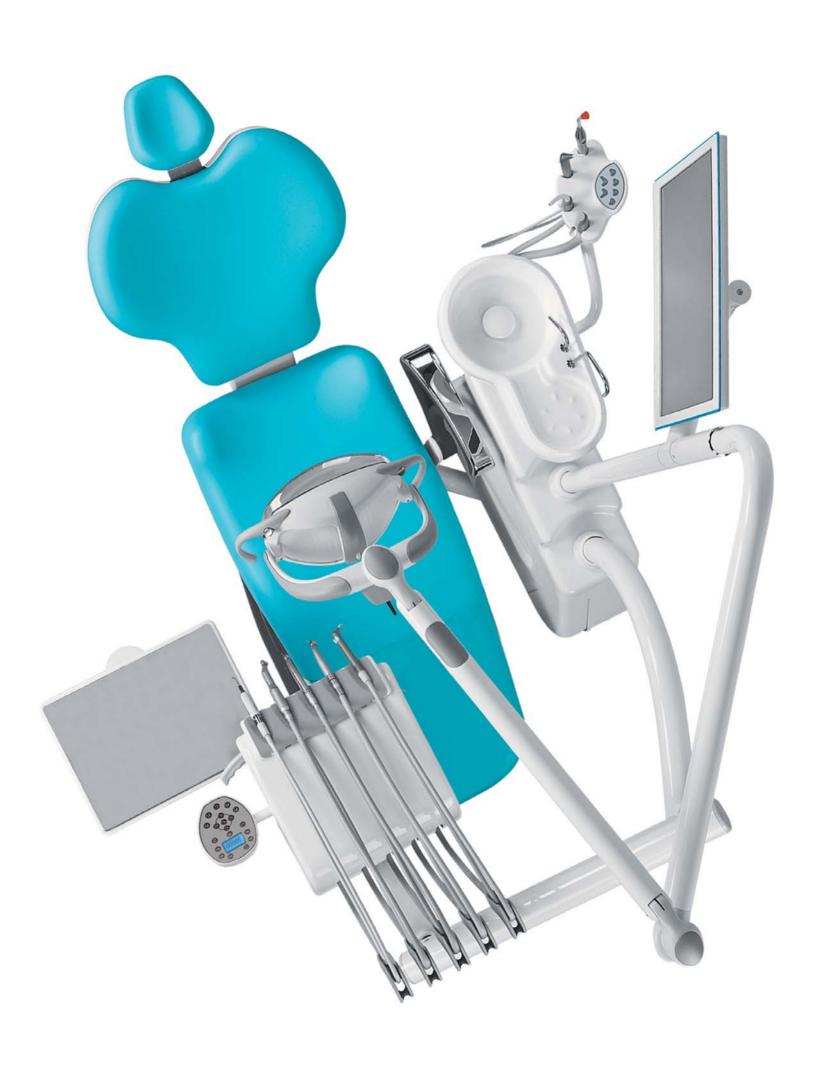

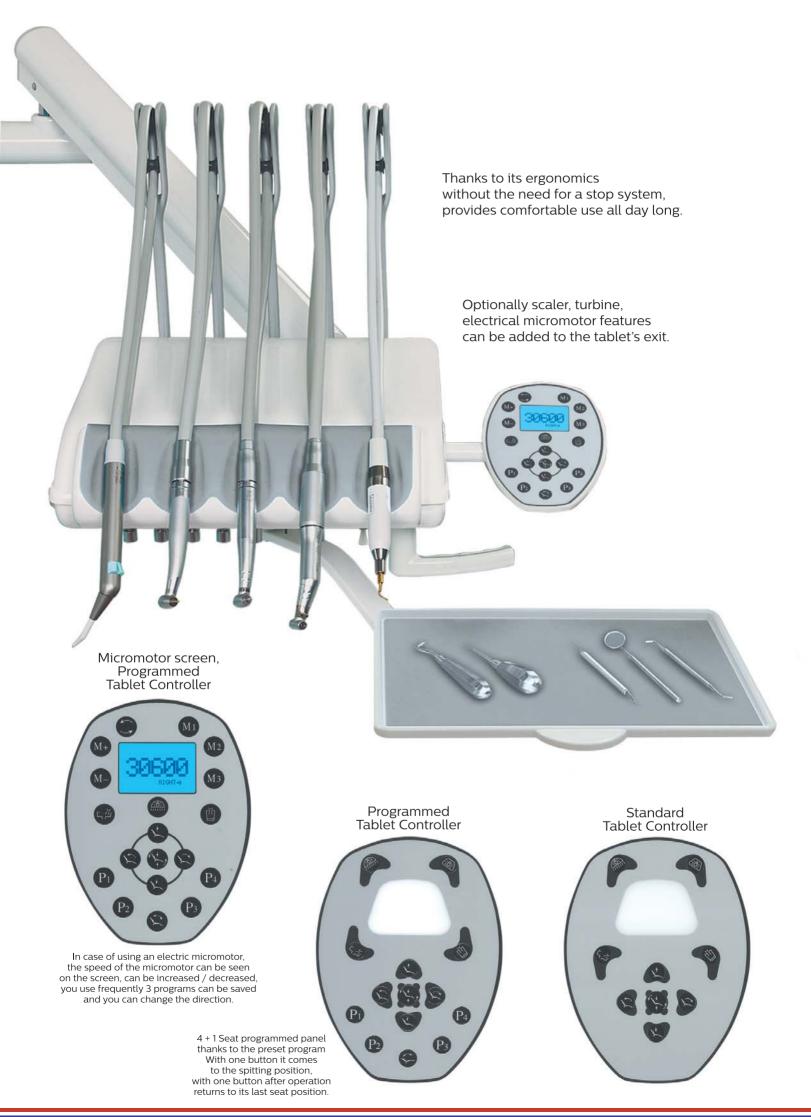



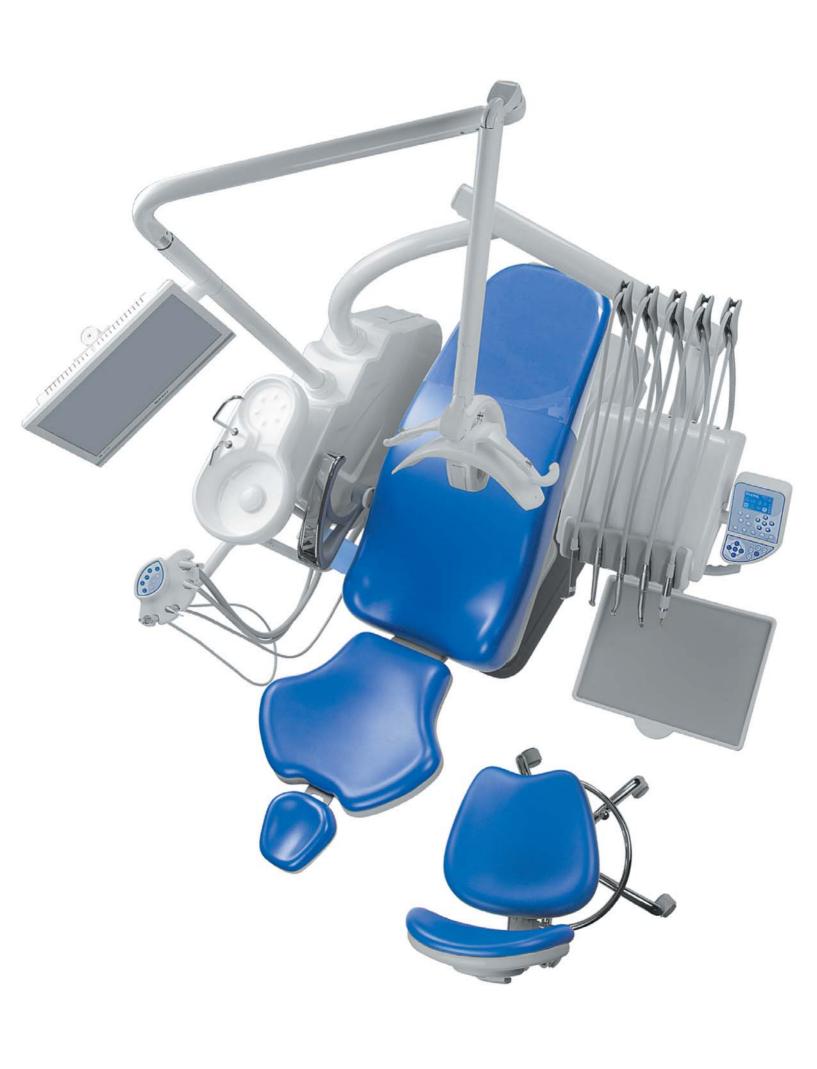







With an easy to use interface, with one button it is possible; cuspidor glass water on /off, selection of water tank or city water or unit disinfection system, working with water or without, peristaltic pump solution speed, chair movement control and chair movement programm.





Optionally in addition to 2 surgical aspiration exit, polymerization, second syringe or other handpieces can be added.





Useful for hygiene and less sensitive against disinfection solutions, it can be removed easily and autoclavable.





Inox cuspidor and glasses pipes make it possible use them easily without having problem such a rust or peeling.





JUST IMAGINE,






WE REALISE YOUR DREAMS





www.tekmil.com.tr

micromotor selection, scaler, control of cuspidor and class with timer

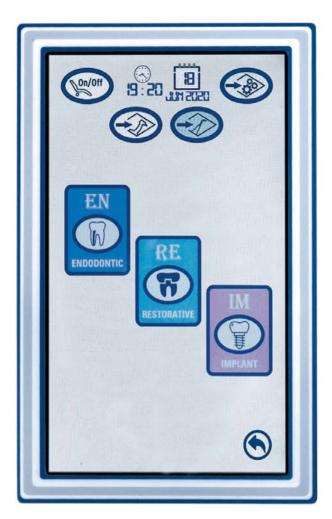



RIGHT OF CHOICE

From micromotor menu, you can select endodontic, restorative and implant micromotor Correcta Unit make it possbile to use double type micromotor at the same unit.

ENDODONTIC SCREEN

A new user friendly screen allows you control of speed, torque, light intensity, liquid level of peristaltic pump, contraangle selection, micromtoor 3 program speed, micromotor 3 program mode (auto reverse, normal mode, auto forward)

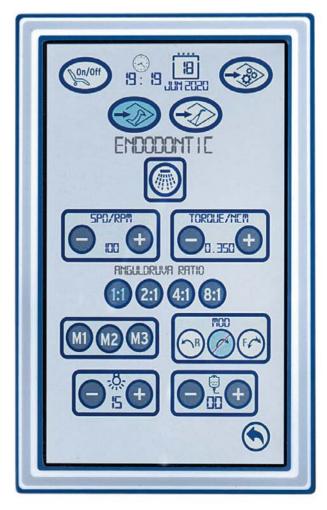

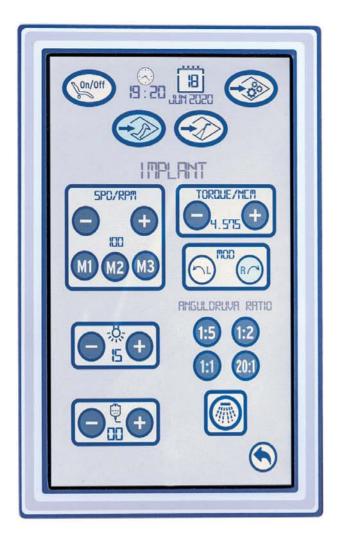

IMPLANT

Touch panel allows you control in addition to speed , also torque , light intensity, micromotor 3 program speed, direction of the micromotor (left and right) liquid level of peristaltic pump, contraangle selection, In addition, with one button it is possible turn on / off water

RESTORATIVE

The micromotor speed can be controlled between 1000 -40.000 rpm, you can select the contraanlge and control, light intensity, liquid level of peristaltic pump, contraangle selection, micromotor 3 program speed, direction of the micromotor (left and right)





SCALER SCREEN

Programmed Endo, Perio, Scaler 3 program allows you to make proper tooth cleaning and with one button it is possible to turn on/off the water





STEP BY STEP IN HYGIENE





Aspiration filters and tubes are easy to remove and autoclavable.





UNIT DISINFECTION SYSTEM

With %3 hydrogen peroxide, this system make it possible the cleaning of the whole way to glasses and tablet's exits.



In case of need of city water, you can just select and continue to your treatment.





Peristaltic Pump Provide sterile irrigation liquid to micromotors for surgical operation.





Magnetical toolhandle place, easy to remove and easy to clean.













Narrow or wide Backrest depends on your working habits



Trandelenburg movement can be added to the chair movement, standart or with programmed version.



Movable **right arm** can be added optionally to have more comfort for your patients.

Control that can be fixed to the base of the unit, gives user possibility to change chair movement without using hands.





Multi axis movement

Freedom in the working position Left right, in every angle; possibility to fix the headrest.





Double Movement Headrest

With double movement dentist can adjust easily most comfortable position for patient. For children adjust down,

for wheelchair patient easy in reverse position.






Multifunction chrome based foot rest dentist stool with automatic brake



The backrest part is movable up/ down and forward / backward.

It is possible to fix in desired position.

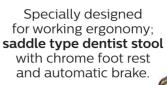





Electrical pneumatical pedal with chair movement control (chip blower system)








Pneumatic pedal Chip Blower system, with water spray control



What is your COLOR of life?
Wide range selection of colors suitable for your working environment



Correcta Unit Dimensions

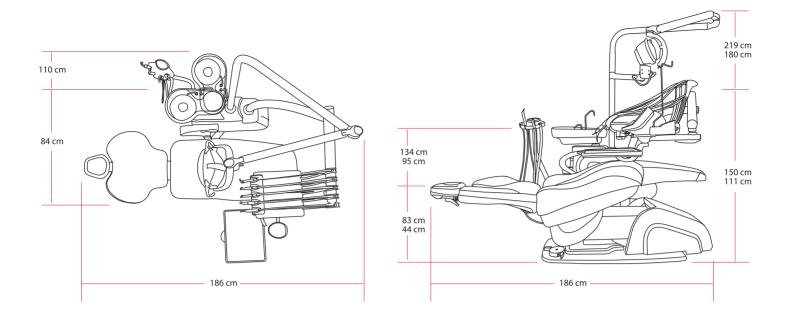









YOUR HEALTH IS IMPORTANT FOR US Therefore we give importance HYGIENE with Correcta Unit.





💡 İnönü Mahallesi Çayır Sokak No:43 Harbiye-Şişli-İstanbul/ TURKEY

+90 212 233 7837

+90 212 231 7277 × satis@tekmil.com.tr

230 4986

+90 545 96367 26



www.tekmil.com.tr